

Canine Haptoglobulin ELISA Kit

Vertrieb:

 L O X O GmbH
 Immunbiologie Biochemie, Produkte und Systeme

 Postfach 11 30
 69215 Dossenheim

 Telefon +49 (0) 62 21 - 86 80 23
 FAX +49 (0) 62 21 - 86 80 255

 E-Mail: info@loxo.de
 Internet: www.loxo.de

Assaypro LLC 30 Triad South Drive St. Charles, MO 63304 T (636) 447-9175 F (636) 447-9475

www.assaypro.com


Hinweis/Note: Der Packungsbeileger dient nur als erste Information. Der relevante Packungsbeileger liegt der Ware bei.

The datasheet is only a first information. The relevant datasheet is included with the product.

For any questions regarding troubleshooting or performing the assay, please contact our support team at support@assaypro.com.

Thank you for choosing Assaypro.

Assay Summary

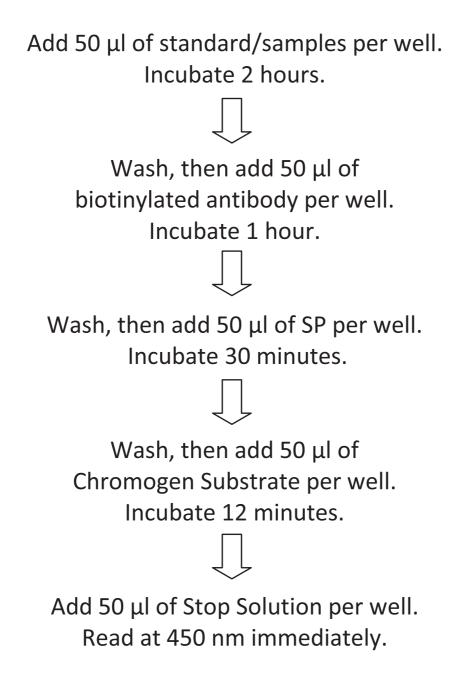

	1	2	3	4	5	6	7	8	9	10	11	12
A												
В												
С												
D												
E												
F												
G												
н												

Assay Template

AssayMax Canine Haptoglobulin ELISA Kit

Catalog No. ECH2003-1 Sample Insert/Reference Only

Introduction

Haptoglobulin (Hpt) is a plasma protein with hemoglobin-binding capacity, and a plasma glycoprotein that forms a stable complex with hemoglobin to aid the recycling of heme iron. It is a well-known marker of hemolysis (1). High haptoglobulin level in plasma is associated with an increased cardiovascular risk in obese men (2), inflammation (3), atherosclerosis (4), and systemic sclerosis (5).

Principle of the Assay

The AssayMax Canine Haptoglobulin ELISA kit is designed for detection of canine haptoglobulin in plasma and cell culture supernatants. This assay employs a quantitative sandwich enzyme immunoassay technique that measures haptoglobulin in less than 4 hours. A polyclonal antibody specific for haptoglobulin has been pre-coated onto a microplate. Haptoglobulin in standards and samples is sandwiched by the immobilized antibody and a biotinylated polyclonal antibody specific for haptoglobulin, which is recognized by a streptavidin-peroxidase conjugate. All unbound material is then washed away and a peroxidase enzyme substrate is added. The color development is stopped and the intensity of the color is measured.

Caution and Warning

- Prepare all reagents (working diluent buffer, wash buffer, standards, biotinylated antibody, and SP conjugate) as instructed, prior to running the assay.
- Prepare all samples prior to running the assay. The dilution factors for the samples are suggested in this protocol. However, the user should determine the optimal dilution factor.
- Spin down the SP conjugate vial and the biotinylated antibody vial before opening and using contents.
- This kit is for research use only.
- The kit should not be used beyond the expiration date.
- The Stop Solution is an acidic solution.

Reagents

- **Canine Haptoglobulin Microplate:** A 96-well polystyrene microplate (12 strips of 8 wells) coated with a polyclonal antibody against canine haptoglobulin.
- **Sealing Tapes:** Each kit contains 3 precut, pressure sensitive sealing tapes that can be cut to fit the format of the individual assay.
- **Canine Haptoglobulin Standard:** Canine Haptoglobulin in a buffered protein base (40 ng, lyophilized).
- **Biotinylated Canine Haptoglobulin Antibody (50x):** A 50-fold concentrated biotinylated polyclonal antibody against canine haptoglobulin (140 µl).
- **MIX Diluent Concentrate (10x)**: A 10-fold concentrated buffered protein base (30 ml).
- Wash Buffer Concentrate (20x): A 20-fold concentrated buffered surfactant (30 ml, 2 bottles).
- Streptavidin-Peroxidase Conjugate (SP Conjugate): A 100-fold concentrate (80 μl).
- **Chromogen Substrate**: A ready-to-use stabilized peroxidase chromogen substrate tetramethylbenzidine (8 ml).
- **Stop Solution**: A 0.5 N hydrochloric acid to stop the chromogen substrate reaction (12 ml).

Storage Condition

- Store components of the kit at 2-8°C or -20°C upon arrival up to the expiration date.
- Store SP Conjugate and biotinylated antibody at -20°C.
- Store Microplate, Diluent Concentrate (10x), Wash Buffer, Stop Solution, and Chromogen Substrate at 2-8°C.
- Unused microplate wells may be returned to the foil pouch with the desiccant packs and resealed. May be stored for up to 1 month in a vacuum desiccator.
- Diluent (1x) may be stored for up to 1 month at 2-8°C.
- Store standard at 2-8°C before reconstituting with diluent and at -20°C after reconstituting with diluent.

Other Supplies Required

- Microplate reader capable of measuring absorbance at 450 nm.
- Pipettes (1-20 μl, 20-200 μl, 200-1000 μl and multiple channel).
- Deionized or distilled reagent grade water.

Sample Collection, Preparation and Storage

- **Cell Culture Supernatants:** Centrifuge cell culture media at 3000 x g for 10 minutes to remove debris. Collect supernatants and assay. Store the remaining samples at -20°C or below. Avoid repeated freeze-thaw cycles
- **Plasma:** Collect plasma using one-tenth volume of 0.1 M sodium citrate as an anticoagulant. Centrifuge samples at 3000 x g for 10 minutes. Dilute samples 1:200000 with MIX Diluent and assay. The undiluted samples can be stored at -20°C or below for up to 3 months. Avoid repeated freeze-thaw cycles. (EDTA or Heparin can also be used as an anticoagulant.)

Reagent Preparation

- Freshly dilute all reagents and bring all reagents to room temperature before use.
- MIX Diluent Concentrate (10x): If crystals have formed in the concentrate, mix gently until the crystals have completely dissolved. Dilute the MIX Diluent Concentrate 1:10 with reagent grade water. Store for up to 1 month at 2-8°C.
- Standard Curve: Reconstitute the 40 ng of Canine Haptoglobulin Standard with 2 ml of MIX Diluent to generate a solution of 20 ng/ml. Allow the standard to sit for 10 minutes with gentle agitation prior to making dilutions. Prepare duplicate or triplicate standard points by serially diluting the standard solution (20 ng/ml) 1:2 with MIX Diluent to produce 10, 5, 2.5, 1.25, 0.625, and 0.313 ng/ml solutions. MIX Diluent serves as the zero standard (0 ng/ml). Any remaining solution should be frozen at < -20°C.

Standard Point	Dilution	[Canine Haptoglobulin] (ng/ml)		
P1	Standard (20 ng/ml)	20.00		
P2	1 part P1 + 1 part MIX Diluent	10.00		
P3	1 part P2 + 1 part MIX Diluent	5.000		
P4	1 part P3 + 1 part MIX Diluent	2.500		
P5	1 part P4 + 1 part MIX Diluent	1.250		
P6	1 part P5 + 1 part MIX Diluent	0.625		
P7	1 part P6 + 1 part MIX Diluent	0.313		
P8	MIX Diluent	0.000		

• **Biotinylated Canine Haptoglobulin Antibody (50x):** Spin down the antibody briefly and dilute the desired amount of the antibody 1:50 with MIX Diluent. Any remaining solution should be frozen at -20°C.

- Wash Buffer Concentrate (20x): If crystals have formed in the concentrate, mix gently until the crystals have completely dissolved. Dilute the Wash Buffer Concentrate 1:20 with reagent grade water.
- **SP Conjugate (100x):** Spin down the SP Conjugate briefly and dilute the desired amount of the conjugate 1:100 with MIX Diluent. Any remaining solution should be frozen at -20°C.

Assay Procedure


- Prepare all reagents, working standards and samples as instructed. Bring all reagents to room temperature before use. The assay is performed at room temperature (20-30°C).
- Remove excess microplate strips from the plate frame and return them immediately to the foil pouch with desiccants inside. Reseal the pouch securely to minimize exposure to water vapor and store in a vacuum desiccator.
- Add 50 μ l of Canine Haptoglobulin Standard or sample per well. Cover wells and incubate for 2 hours. Start the timer after the last sample addition.
- Wash five times with 200 µl of Wash Buffer manually. Invert the plate each time and decant the contents; hit 4-5 times on absorbent material to completely remove the liquid. If using a machine, wash six times with 300 µl of Wash Buffer and then invert the plate, decanting the contents; hit 4-5 times on absorbent material to completely remove the liquid.
- Add 50 μl of Biotinylated Canine Haptoglobulin Antibody to each well and incubate for 1 hour.
- Wash the microplate as described above.
- Add 50 μl of Streptavidin-Peroxidase Conjugate per well and incubate for 30 minutes. Turn on the microplate reader and set up the program in advance.
- Wash the microplate as described above.
- Add 50 µl of Chromogen Substrate per well and incubate for about 12 minutes or till the optimal blue color density develops. Gently tap the plate to ensure thorough mixing and break the bubbles in the well with pipette tip.
- Add 50 μl of Stop Solution to each well. The color will change from blue to yellow.
- Read the absorbance on a microplate reader at a wavelength of 450 nm immediately. If wavelength correction is available, subtract readings at 570 nm from those at 450 nm to correct optical imperfections. Otherwise, read the plate at 450 nm only. Please note that some unstable black particles may be generated at high concentration points after stopping the reaction for about 10 minutes, which will reduce the readings.

Data Analysis

- Calculate the mean value of the duplicate or triplicate readings for each standard and sample.
- To generate a standard curve, plot the graph using the standard concentrations on the x-axis and the corresponding mean 450 nm absorbance on the y-axis. The best-fit line can be determined by regression analysis using log-log or four-parameter logistic curve-fit.
- Determine the unknown sample concentration from the Standard Curve and multiply the value by the dilution factor.

Standard Curve

• The curve is provided for illustration only. A standard curve should be generated each time the assay is performed.

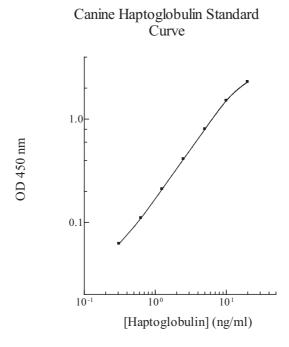

Performance Characteristics

- The minimum detectable dose of canine haptoglobulin is typically ~ 0.3 ng/ml.
- Intra-assay and inter-assay coefficients of variation were 5.1% and 7.1 % respectively.

Linearity

	Average Percentage of Expected Value
Sample Dilution	Plasma
1:100000	83%
1:200000	97%
1:400000	104%

Recovery

Standard Added Value	1 – 10 ng/ml
Recovery %	83-118 %
Average Recovery %	96 %

Cross-Reactivity

Species	% Cross Reactivity		
Human	<1%		
Bovine	None		
Monkey	<5%		
Mouse	None		
Rat	None		
Swine	None		
Canine	100%		

• 10% FBS in culture media will not affect the assay.

References

- (1) Van Vlierberghe H *et al* (2004) *Clin Chim Acta*. 345(1-2): 35-42
- (2) Engstrom G et al. (2004) Arterioscler Thromb Vasc Biol. 24(8): 1498-502
- (3) Rocha-Pereira P et al. (2004) Br J Dermatol. 150(5): 917-28
- (4) Matuszek MA et al. (2003) Atherosclerosis 168 (2): 389-96
- (5) Kucharz EJ et al. (2000) Clin Rheumatol 19(2):165-6

Version 1.4

Related Products

- EH1003-1 AssayMax Human Haptoglobulin ELISA Kit (Plasma and Serum samples)
- EH2003-1 AssayMax Human Haptoglobulin ELISA Kit (Urine, Milk, Saliva, Cell Culture samples)
- ERH1003-1 AssayMax Rat Haptoglobulin ELISA Kit (Plasma and Serum samples)
- ERH2003-1 AssayMax Rat Haptoglobulin ELISA Kit (Urine and Cell Culture samples)
- EBH2003-1 AssayMax Bovine Haptoglobulin ELISA Kit (Urine and Cell Culture samples)
- EPH2003-1 AssayMax Swine Haptoglobulin ELISA Kit (Plasma, Serum, Urine, and Cell Culture samples)